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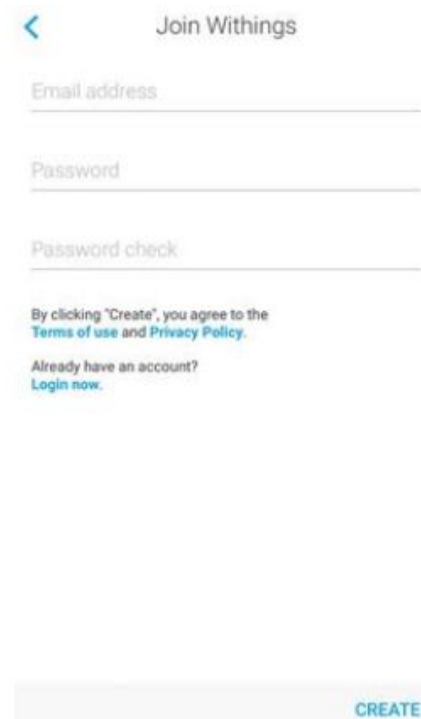
The screenshot shows a mobile app interface for joining Withings. At the top left is a back arrow, and at the top right is the title 'Join Withings'. Below the title are three input fields: 'Email address', 'Password', and 'Password check'. Underneath these fields is a line of text: 'By clicking "Create", you agree to the [Terms of use and Privacy Policy](#).' Below that is another line: 'Already have an account? [Login now](#).' At the bottom right of the form is a blue button labeled 'CREATE'.

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
The screenshot shows a mobile app interface for the 'Profile' registration step. At the top center is the title 'Profile'. Below the title are four input fields: 'First name', 'Last name', 'Sex', and 'Date of birth'. The 'Sex' field has two circular icons representing male and female. At the bottom right of the form is a blue button labeled 'NEXT'.

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The screenshot shows a mobile application interface for profile registration. At the top left is a back arrow icon, and at the top center is the title "Profile". Below the title are two input fields: "Height" and "Weight", each with a horizontal line underneath. At the bottom right of the screen is a light gray button with the text "NEXT" in blue.

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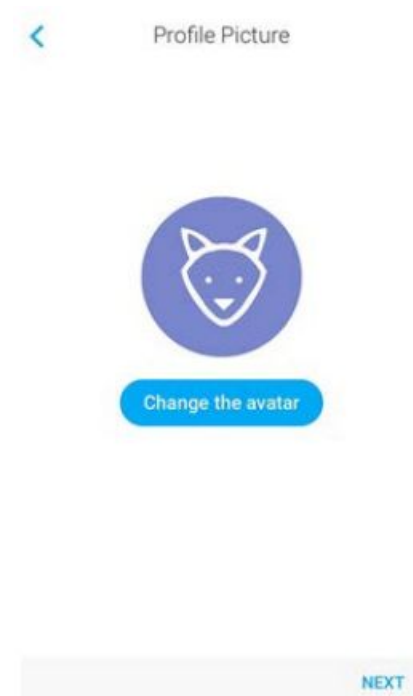The screenshot shows a mobile application interface for selecting a profile picture. At the top left is a back arrow icon, and at the top center is the title "Profile Picture". In the center of the screen is a circular blue icon containing a white outline of a cat's face. Below this icon is a blue button with the text "Change the avatar" in white. At the bottom right of the screen is a light gray button with the text "NEXT" in blue.

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
Cancel Installation

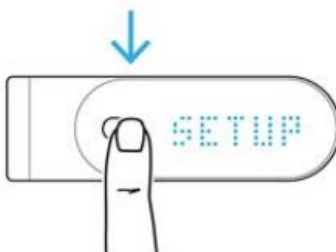Press and hold the button until the display shows "SETUP".

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